

Investigator™ 12 Green

Single Use Gamma Irradiated Hollow Fiber Cartridges. Easy and Reliable

The Investigator™ Green crossflow hollow fiber membrane cartridges are designed for pilot scale and small scale production and is the perfect tool for process development. With a surface area of 0.13 m², the Investigator™ 12 Green can accommodate batch volumes that range from 1L to 6L with an average permeate flow rate of 3.9 to 7.8 L/hr for ultrafiltration applications and 2.6 to 5.2 L/hr for microfiltration applications.


The Investigator™ Green cartridges are gamma irradiated and ready to use, and utilize WaterSep's well accepted antifouling, low binding, modified polyethersulfone membrane (m-PES) which provides a process flux and product recovery that typically exceeds other crossflow devices.

The Investigator™ 12 Green hollow fiber cartridges are offered with the same path length as other high performance WaterSep HF cartridges, which makes scale up/scale down easy and predictable.

Investigator™ Green cartridges are also offered as part of fully integrated open architecture single use crossflow assemblies available with pressure sensors, flow meters, pumpheads, aseptic connectors or bags.

The Investigator™ 12 Green cartridges offer a true single use crossflow system that provides:

- Significant cost savings over traditional cassette formats.
- Gamma irradiated - ready to use - no pre-rinse required.
- No humectants - 80x less extractables, see [Extractables Flush Analysis](#).
- No need for time-consuming cleaning validation.
- Reduced risk of cross contamination.
- Consistent membrane performance batch to batch.
- Self-contained – no assembly – no expensive hardware – no membrane installation.
- Every HF cartridge is integrity tested and has an individual lot number for easy traceability.

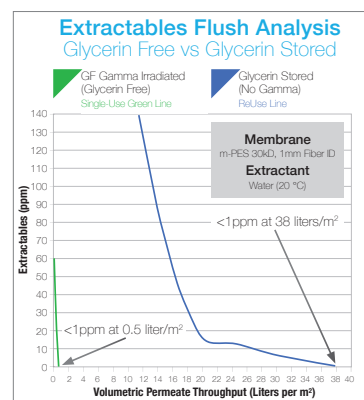

TYPICAL APPLICATIONS

- Clarification of **mammalian/CHO cell** cultures, providing maximum protein recovery.
- **Cell-harvest**. Excellent results have been achieved with both *E. coli* whole cells and *E. coli* lysates, as well as other microbial process streams.
- Concentration and purification of **vaccines**.
- Concentration and diafiltration of **gene therapy** products.
- Concentration/diafiltration of **monoclonal antibodies, recombinant proteins, biological macromolecules** and **peptides**.

Investigator™ 12 Green

ORDERING INFORMATION

Place orders online at www.pure-gmbh.com, or email purchase orders to support@pure-gmbh.com, sales order confirmation to follow.

Investigator™ 12 Green HF Cartridge

0.5 mm Fiber ID

SU XXX 05INV12 S3, 3/pkg

SU XXX 05INV12 S6, 6/pkg

1.0 mm Fiber ID

SU XXX 10INV12 S3, 3/pkg

SU XXX 10INV12 S6, 6/pkg

2.0 mm Fiber ID*

SU XXX 20INV12 S3, 3/pkg

SU XXX 20INV12 S6, 6/pkg

*Available in 10K, 50K, 300K-750K and 0.1µm-0.65µm

Replace XXX with:

003 for 3K	100 for 100K	910 for 0.1µm
005 for 5K	300 for 300K	920 for 0.2µm
010 for 10K	500 for 500K	945 for 0.45µm
030 for 30K	750 for 750K	965 for 0.65µm
050 for 50K		

SPECIFICATIONS

Cartridge Dimensions

1.32" (3.34 cm) x 12.01" (30.50 cm)

Membrane Surface Area

1.4 ft² (0.13 m²)

Molecular Weight Cut-off

3K, 5K, 10K, 30K, 50K, 100K, 300K, 500K, 750K

Membrane Pore Size

0.1µm, 0.2µm, 0.45µm and 0.65µm

Fiber ID

0.5 mm, 1.0 mm, 2.0 mm

MATERIALS OF CONSTRUCTION

The cartridge as assembled is USP Class VI compliant.

Membrane

Modified Polyethersulfone (m-PES)

Housing

White Polysulfone

Encapsulant

Medical Grade Epoxy


CONNECTIONS

Feed/Retentate

1.5" TC Standard Sanitary Fitting

Permeate

0.5" TC Standard Sanitary Fitting

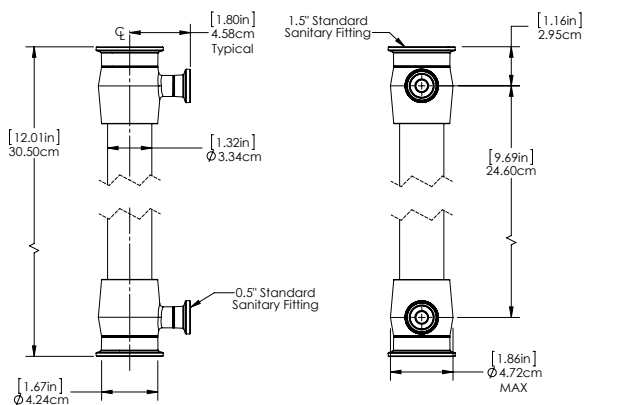

Technical Support: support@pure-gmbh.com
or call +49 6431 49612 40.